
SOGO

Human Technology

manual de instrucciones / instructions for use
mode d'emploi / manual de instruções

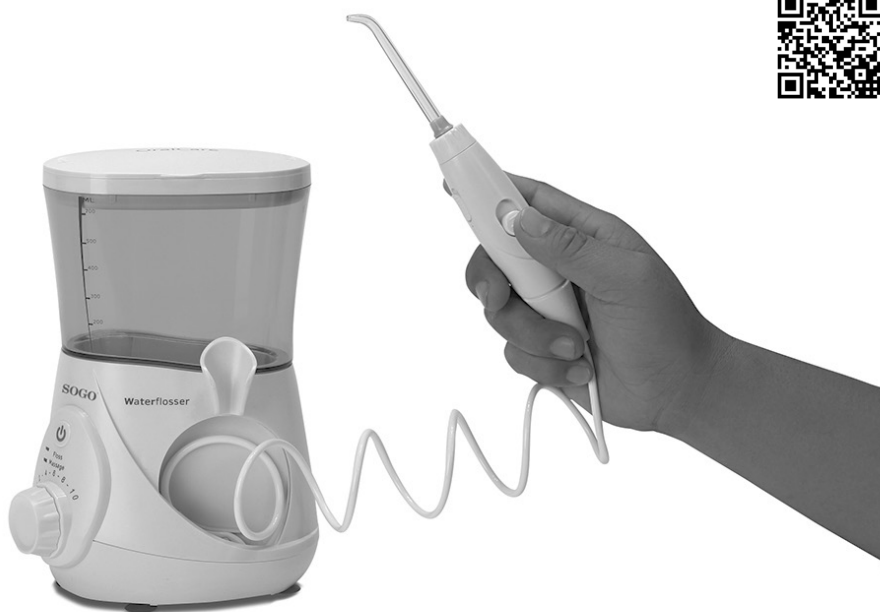Irrigador Dental y Bucal Water Flosser Flosser de l'eau Irrigadores Bucais

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Potencia de Consumo: 50W

Instructions for use

- Welcome! By selecting the SOGO Water Flosser you are on your way to have better gum health, and a healthier smile!
- With the combination of water pressure and pulsations, the SOGO Water Flosser removes harmful bacteria deep between the teeth and below the gumline where brushing alone won't reach.

Traditional Oral Clean

- Traditional brushing can't easily reach hidden food debris and harmful bacteria.
- String Flossing has limited access below the gumline

Important Safeguards

- When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

This manual can also be downloaded from our website www.sogo.es

Read all instructions carefully before using.

Danger:

To reduce the risk of electrocution:

- Always unplug product after using.
- Do not handle plug with wet hands.
- Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- Do not use while bathing.
- Do not immerse in water or other liquid.


- Do not reach for a product that has fallen into any liquid. Unplug immediately.
- Check the power supply cord for damage before the first use and during the life of the product.
- Only use attachments, including jet tips Tooth Pocket tips, tongue cleaner tips and orthodontic tips that are recommended by us.
- Do not remove the Tooth Pocket flexible tip from the clear shaft. If it separates for any reason, discard both the tip and shaft and replace with a new Tooth Pocket tip.
- Fill reservoir with water or other recommended solutions only.

Warning:


- To reduce the risk of burns, electric shock, fire, serious injury or damage to tissue:
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Return the product to the local distributor.
- Closely instruct and supervise children and individuals with special needs in the proper use of this product.
- Do not direct water under the tongue, into the ear, nose, or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas. See “operating instructions” for correct usage.
- Use this product only as indicated in these

instructions or as recommended by your dentist.

- Do not use this product for more than 5 minutes in each two-hour period.
- Do not drop or insert any foreign object into any opening or hose.
- Keep electrical cord away from heated surfaces.
- Do not operate where oxygen or aerosol sprays are being used.
- Do not use outdoors.
- Do not use while wearing any oral jewelry. Remove prior to use.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- Do not use this product for more than 5 minutes in each two-hour period.
- Children should be supervised to ensure that they do not play with the appliance.
- If the power cord is damaged, it must be replaced by the manufacturer or its service agent or a similarly qualified person in order to avoid a hazard

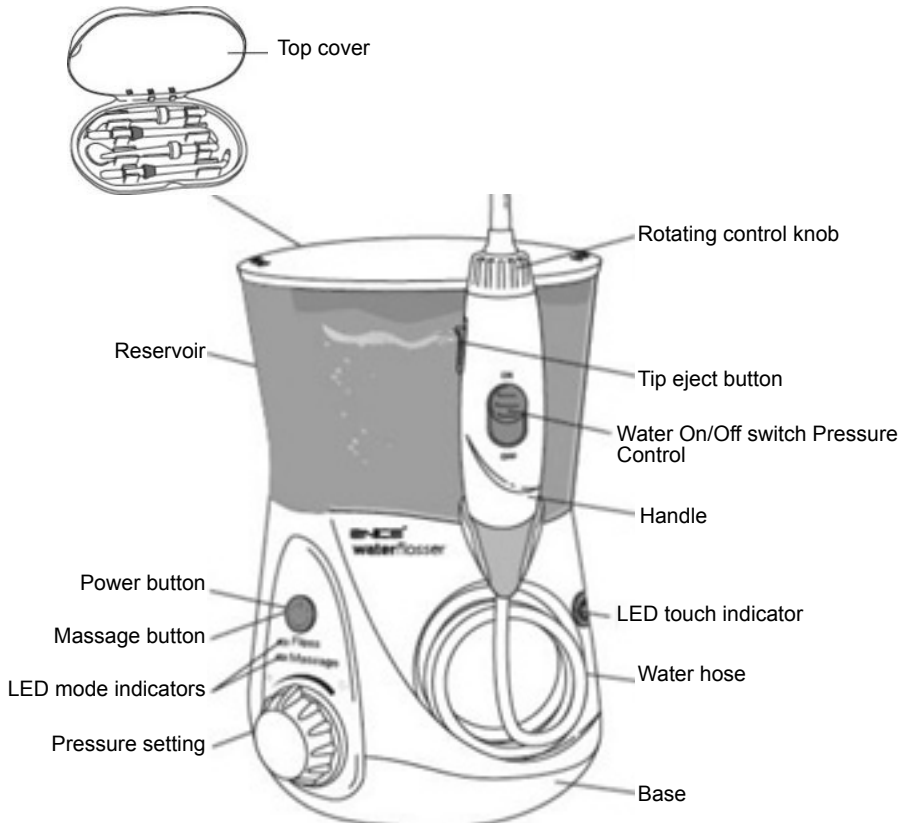
Jet Tip

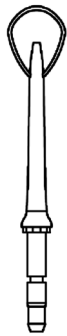
Tongue
Cleaner Tip

Pik Pocket
Tip

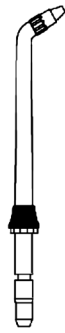Orthodontic
Tip

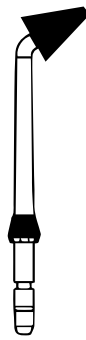
Nasal

OPERATING INSTRUCTIONS


- **Connection**
Plug the cord into an electrical outlet. If the outlet is controlled by a wall switch, be sure it is turned on.
- **Preparing the Reservoir**
Remove reservoir and cover from base. Fill the reservoir with lukewarm water. Place the reservoir over the base and press down firmly.


- **Inserting and removing Tips**
Insert tip into the center of the knob at the top of the dental jet handle. The colored ring will touch the knob if the tip is correctly locked into place. To remove tip from the handle, press the jet tip eject button and pull tip from handle

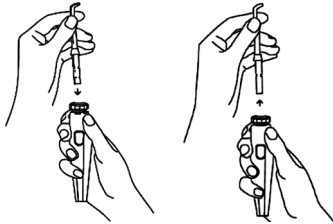
- **Adjusting the Pressure Setting**
Turn the pressure control dial on the base of the unit to the lowest setting for first-time use. Gradually increase pressure over time to the setting you prefer-or as instructed by your dental professional.




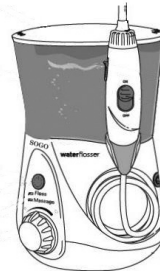
- **Placement in Mouth**
Lean low over sink and place the jet tip in mouth. Aim the tip toward teeth. Turn unit on: press the power button for 3 seconds for flossing mode, press again for pulsing massage mode; with Blue LED indicator represents flossing. Green represents pulsing.
Set pressure at desired setting.

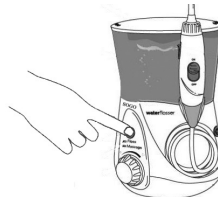Direct the stream by turning the knob at the top of the handle.

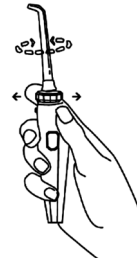- **Recommended Technique**
Direct the jet stream at 90-degree angle to your gumline. Slightly close lips to avoid splashing but allow water to flow freely from mouth into the sink. For best results, start in the molar area (back teeth) working toward the front teeth. Glide jet tip along gumline and pause briefly between teeth. Continue process until all areas around and between the teeth have been cleaned.

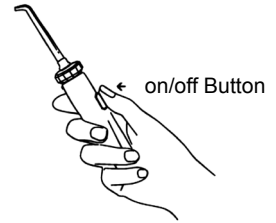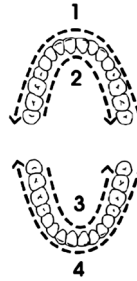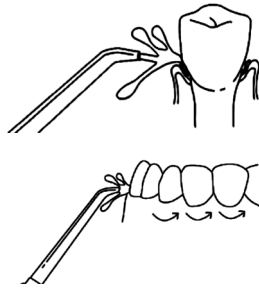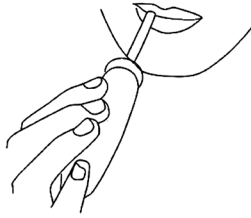On/Off Control

You can stop the flow any time by pressing the on/off button on the handle

When Finished

Turn the unit off. Empty any liquid left in the reservoir. Unplug the unit.


• **Hydro-Pulse Massage Mode**

Hydro-Pulse Massage offers a unique water pulsation pattern designed to enhance gum stimulation. For best results, use the Floss Mode as instructed. After using Floss Mode for 1 minute, switch the unit to Massage Mode, adjust pressure as desired, and aim the water stream at your gums and massage throughout your mouth.

• **Special Tips**

The Tooth Pocket Tip delivers water or anti-bacterial solutions deep into periodontal pockets.


To use the Tooth Pocket Tip, lean low over your sink and place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gum-line, into the pocket. Turn the unit on to the lowest pressure setting and clean by gently following gum-line, inserting tip in between teeth. Continue tracing along the gum-line until you finish cleaning.

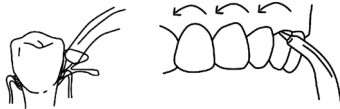• **Orthodontic Tip**

The special Orthodontic Tip is uniquely designed to simultaneously brush and rinse hard-to-reach areas around braces and other dental work.

To use the orthodontic tip, lean low over the sink and use the brush bristles to clean around dental work. Water will flow through the brush to rinse away any debris that is removed.

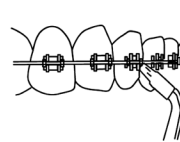Tongue Cleaner Tip

Lean over the sink. Place the tongue cleaner on your tongue and start the unit-use the lowest pressure setting, increase as comfortable. Gently place the tongue cleaner in the center/midline of your tongue about half way back and with light pressure pull forward.

• **Nasal**

Rise away dust, dirt, allergens and pollen. Clears congestion & help prevent sinus infections

USE OF THE TIPS

Type	Used For	Advantage
	Classic Jet Tip General Use	More effective than string floss for healthier gums
	Nasal tip Rinse away dust, dirt, allergens and pollen.	Clears congestion & help prevent sinus infections
	Orthodontic Tip For Braces	3X as effective as string floss for removing plaque around Orthodontic applications
	Pik Pocket Tip Periodontal Pockets - Furcations	Gently delivers medicaments and rinses deep below the gumline
	Tongue Cleaner Fresh Breath	Removes bacteria from the tongue to freshen breath

PROBLEM	CAUSE	CORRECTIVE ACTION
All Tips Leakage from the handle knob.	Tip not fully engaged.	Remove top and reinstall
Tooth Pocket Tip Tear in soft tip.	Tip is worn	Replace with new tip
Water Flosser Base Hose to the handle leaks.	Damage to the hose.	Can be replaced by the technician
Reservoir leaks	Missing reservoir valve or reservoir valve is upside down	Replace reservoir valve Turn valve right-side up
Inadequate pressure	Reservoir not positioned properly	Re-position reservoir. Press down firmly on the reservoir
Unit doesn't start	Motor not running	Be sure outlet is functioning; test with another appliance (check electrical reset button and/or wall switch).

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DECLARATION OF CONFORMITY

This device complies with Electromagnetic Compatibility Directive 2014/30/EU, Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment and Directive 2009/125/EC on the eco-design requirements applicable to energy-related products.

DECLARATION OF CONFORMITY

This device complies with Electromagnetic Compatibility Directive 2014/30/EU, Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment and Directive 2009/125/EC on the eco-design requirements applicable to energy-related products.

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This symbol means that in case you wish to dispose of the product once its working life has ended, take it to an authorised waste agent for the selective collection of waste electrical and electronic equipment (WEEE)

Imported by:

Sanysan Appliances S.L, NIF: B98753056,
C/ Barcas 2, 2
46002 Valencia, Spain
Product made in P.R.C. – Designed by
SOGO based on the European Quality
Standards
Customer Service: www.sogosat.com /
sogosat@sogosat.com / 0034 902 222 161

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