SOGO

Human Technology

manual de instrucciones / instructions for use
mode d'emploi / manual de instruções
gebrauchsanweisung / manuale d'istruzioni

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IMPORTANT SAFEGUARDS

This instruction manual can also be downloaded from our web page www.sogo.es

Read the operating instructions carefully before putting the appliance into operation and keep the instructions including the warranty, the receipt and, if possible the box with the internal packing. If you give this device to other people, please also pass on the operating instruction. It is essential to comply with these instructions in order to avoid accidents and prevent damage to the machine:

- The device is designed exclusively for private and not for commercial use.
- Only use the device for the intended purpose.
- Do not use a damaged device.
- Keep away from heat, direct sunlight, damp (never immerse in liquids) and sharp edges.
- Do not try to repair the appliance on

your own. Always contact an authorized technician.

- In order to ensure your children's safety, please keep all packaging (plastic bags, boxes, polystyrene etc.) out of their reach.
- Children should be supervised at all times in order to ensure that they do not play with device.
- If the surface of the scale is moist or covered with water droplets, please wipe the surface dry as there is a danger of slipping.
- Handle the scale with care, as it is a precise measuring instrument. Do not drop it or jump on it.
- Do not load the appliance with more than 180kg to avoid damaging the device.
- Always store it only in a horizontal position.
- The body fat analysis scale is not suitable for the following persons:

- Children under 10 years and adults over 80 years

- Persons with symptoms of fever, edema, or osteoporosis

- Persons undergoing treatment for

dialysis

- Persons who use a pacemaker
- Persons who take cardiovascular medication
- Pregnant women
- Sports persons who perform over 5 hours of intensive sporting activities each day
- Person with a heartbeat lower than 60 beats / minute

HANDLING BATTERIES

- Do not expose the battery to high temperatures or direct sunlight. Never throw batteries into the fire. There is a danger of explosion!
- Keep batteries away from children and make sure they don't play with it!
- Do not open the batteries by force.
- Avoid contact with metallic objects (Rings, nails, screws etc.) there is a danger of short-circuiting!
- As a result of a short-circuit batteries may heat up considerably or even catch fire. This may result in burns.
- For your safety the battery poles should be covered with adhesive strips during transport.
- If a battery leaks, make sure your eyes doesn't come in direct contact with the acid. If it does, then wash your eyes thoroughly and rinsed out with clear water, later consult a doctor.
- Different battery types like new and used batteries must not be used together.
- Batteries are not to be disposed of together with domestic wastes. Dispose of the used batteries at the relevant collection bin or at your dealer's.


INSERTING THE BATTERIES

If the supplied batteries have been protected for transport with plastic foil,

1. Remove the foil before inserting the batteries into the appliance.
2. Open the battery compartment on the back.
3. Insert 1pcs Lithium battery; please ensure that the polarity is correct (see the Floor of the battery compartment)!
4. Close the battery compartment again.
5. If the device is not used for a long period of time, please take the batteries out in order to prevent leakage of battery acid.
6. If the message "LO" appears on the display, replace the batteries.

LOCATION OF CONTROLS

- 1 Contact surfaces
- 2 Display
- 3 DOWN-button ∨
- 4 UP-button ∧
- 5 SET-button
- 6 Unit

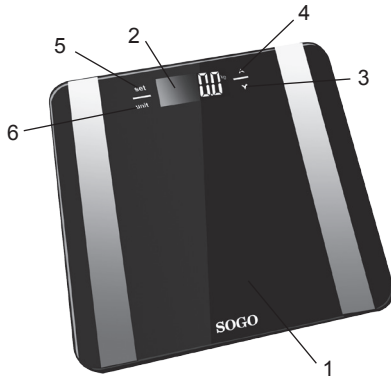INFORMATION ON MEASURING BODY FAT

Do not use the scales if you use a pacemaker or another medical implant with electronic components. The function of the implant may be impaired the measurement.

The device can only function correctly if you stand on the contact surface with bare and dry feet!

The main principle of this device is based on the measurement of the electric impedance within the human body. In order to determine the actual body data, the device reads a weak, unnoticeable electric signal through your body.

This signal is generally safe! However, please read safety instructions!

The method of measurement is referred to as "bioelectric impedance analysis"(BIA), a factor which is based on the relationship between body weight and associated body water percentages as well as on other biological data (age, gender, body weight).

There are heavy persons with a well-balanced percentage between fat and muscle mass, and there are seemingly slim persons with significantly high body fat values. Fat is a factor that puts serious strain on the body, in particular on the cardio-vascular system. Therefore, controlling both values (fat and weight) is a prerequisite for health and fitness.

Display symbol	Meaning
Kg/lb/	Body weight
P0-9	User memory space
	Male
	Female
CM(FT; In)	Body height
Age	Age
FAT	Body fat display
TBW	Body water display
%	Percentage-body fat/ body water/ weight

TECHNICAL SPECIFICATIONS

Power supply: Batteries 2 x CR2032
 Net weight of device: 1.65 kg
 Memory capacity: 10
 Mass alignment: from 75 to 225 cm(2'5"-7'4 1/2"
 "Ft:in)
 Age: from 10 to 80 years
 Minimum weight: from 2.5 kg(5.5lb/0:4st)
 Maximum weight: 150 kg(330.7lb/23:6st)
 Weight measurement unit: 100g increments
 Body fat measuring unit: 0.1%
 Body fat measuring range: 4.0%-50.0%
 Body water measuring range: 35.0%-85.0%
 Technical specifications are subject to changes without prior notice!

BODY FAT

Body fat is an important component of the body. From protecting the joints as a vitamin repository to regulating the body temperature, fat plays a vital role in the body. Thus the goal is not radically reduce the proportion of body fat, but rather to achieve and maintain a healthy balance between muscle and fat.

GUIDELINE VALUES FOR BODY FAT AND BODY WATER CONTENT AS A PERCENTAGE (APPROXIMATE FIGURES)

Age (Years)	Women		Men		
	Body fat	Body water	Body fat	Body water	
< = 30	4.0 - 20.5	66.0 - 54.7	4.0 - 15.5	66.0 - 58.1	Underweight
	20.6 - 25.0	54.6 - 51.6	15.6 - 20.0	58.0 - 55.0	Normal
	25.1 - 30.5	51.5 - 47.8	20.1 - 24.5	54.9 - 51.9	Overweight
	30.6 - 50.0	47.7 - 34.4	24.6 - 50.0	51.8 - 34.4	Obese
> 30	4.0 - 25.0	66.0 - 51.6	4.0 - 19.5	66.0 - 55.4	Underweight
	25.1 - 30.0	51.5 - 48.1	19.6 - 24.0	55.3 - 52.3	Normal
	30.1 - 35.0	48.0 - 44.7	24.1 - 28.5	52.2 - 49.2	Overweight
	35.1 - 50.0	44.6 - 34.4	28.6 - 50.0	49.1 - 34.4	Obese

USE

The device has a very sensitive sensor. If you want to adjust the settings, simply touch the buttons gently with your finger. If the device does not respond to touch, clean the glass surface with a soft, dry cloth. The scales are factory set with the units "kg" and "lb". To change the settings to "lb" or "kg" and "Ft:in", please touch the "unit " button

flash. Confirm the setting by pressing SET or use the UP/DOWN button to choose the symbol for your gender. Then press SET to go to the next setting.

- The height setting will flash. Select your height by pressing UP or DOWN button. Holding down one of the buttons will accelerate the process. Confirm your entry by pressing SET.
- The setting AGE will flash. Select your age by pressing UP or DOWN button. Holding down one of the buttons will accelerate the process. Confirm your entry by pressing SET.

WEIGHING FUNCTION

- Place the scale on a hard, flat surface (avoid carpets). An uneven floor contributes to inaccurate measurements.
- Stand on the scales with both feet. The device switches on automatically.
- Distribute your weight equally and wait for the scale to determine your weight.
- The digits on the display will flash your exact weight.
- In order to extend the battery life, the device will automatically turn off in 10 seconds after you have stepped off the scale.

NOTE

- The device can only function correctly if you stand on the device with bare, dry feet!
- If necessary wipe the contact surface improve results.
- Due to the significant fluctuation of body fat, it is difficult to attain a standard value. On average the optimal testing time is between 6 p.m. and 8 p.m. since the body fat remains very consistent during this period of time.
- Press the SET button. The preset stored position data appears on the display. Use the UP/DOWN button to select the stored position at which your data is saved.
- Wait for the USER display and 0.0kg(0.0lb/st) to stop flashing.
- Stand on the contact surfaces(1) of the scales with both feet.
- Firstly your weight will flash on the display.
- After a few seconds, the following results will be shown:
 -the body fat content in%(FAT symbol)
 -the body water content in %(TBW symbol)
 -the results analysis in the form of a bar

SAVING PERSONAL DATA

- The scale has default values that you can change, as described below.
- Carry out the following steps in sequence and within a few seconds. Otherwise, the device will automatically cancel input mode.
- Press the button SET. The memory space will flash on the display (USER0). Confirm the setting by pressing SET or select another memory space by pressing the UP/DOWN button. Then press SET to go to the next setting.
- The symbol for "male" or "female" will

DISPLAY ERROR MESSAGE

NOTE

If an error message appears on the display, wait until the device is turned off. Turn it on again to restart the operation. If the error display has not disappeared, take the batteries out of the battery compartment for 10 seconds.

Display symbol	Meaning
8888	Self test (when witching on) to check that the numeric display functions properly
----	The Scale stands incorrectly
Lo	Battery is low
ErrL	The percentage of body fat content lies above or below the limits of the measuring scale
Err	The weight lies above the allowed maximum weight

CLEANING

Do not immerse the scale in water: the electronic components can be damaged.
 Do not use a wire brush or any abrasive items to clean the surface.
 Do not use any acidic or abrasive detergents to clean the surface
 Clean the device with a slightly moist cloth. Make sure that no humidity or waters seeps into the device.
 Use a commercially sold detergent or soap for cleaning.
 Dry the device well after cleaning with a dry and soft cloth.

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DECLARATION OF CONFORMITY

This device complies with Electromagnetic Compatibility Directive 2014/30/EU, Directive 2015/863/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment and Directive 2009/125/EC on the eco-design requirements applicable to energy-related products.

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This symbol means that in case you wish to dispose of the product once its working life has ended, take it to an authorised waste agent for the selective collection of waste electrical and electronic equipment (WEEE)

Imported by:

Sanysan Appliances S.L, NIF: B98753056,
C/ Barcas 2, 2
46002 Valencia, Spain
Product made in P.R.C. – Designed by
SOGO based on the European Quality
Standards
Customer Service: www.sogosat.com /
sogosat@sogosat.com / 0034 902 222 161

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